

model **6550**
**WHEELCHAIR
SCALES**



UNIQUE FOLD-UP DESIGN



Detecto
Model 6550
Portable
Wheelchair Scale
with Column

- 1,000 lb / 450 kg capacity
- Lightweight portable fold-up scale
- 2-way low-profile platform for easy patient access



EFFICIENT MEASUREMENT



The 6550's digital weight indicator is the clinical industry's most advanced, featuring EMR/EHR-ready serial, USB and optional Ethernet interface, BMI calculation, yellow "quick keys" for basic weighing, and a multi-line LCD display. The scale is powered by 6 C batteries (not included) or optional AC adapter.



USA-Made Quality

- Up to 99 tares (optional) to store wheelchair weights for efficient patient weighing.
- Hold/Release key to hold weight after leaving the scale for better patient care.
- Interface with a PC or printer via serial, USB or optional Ethernet ports for recordkeeping.
- Patient ID Barcode Scanning: Send 9-digit patient ID to indicator for EMR/EHR records.

One of the most popular medical scales in the industry, this ingeniously-designed USA-made scale offers a wide range of versatility. It can be used with wheelchairs, straight-legged chairs or as a stand-on scale. The 6550's 30 x 32 in / 76 x 81 cm platform (flat portion) features integral sloped ramps allowing access from either side. This lightweight digital fold-up scale is made for heavyweight tasks in patient care.

The 6550's column folds down and latches securely for easy transport in four directions.

Optimum patient weighing "on-the-move"



Integral dual-access ramps allow maximum convenience for loading and unloading patients from either side of the scale.

4-WAY EASY TRANSPORT

The 6550 scale's column folds down and latches securely for easy transport in four directions made possible by the dual-direction wheels. Detecto's 6550 vertical fold-up feature also equals minimal floor space for storage in an upright position when not in use.



Ideal for use in multiple areas



Detecto's 6550 scale features pushbutton and keypad tare for removing the wheelchair weight to display the actual patient weight.



model **6550** WHEELCHAIR SCALES

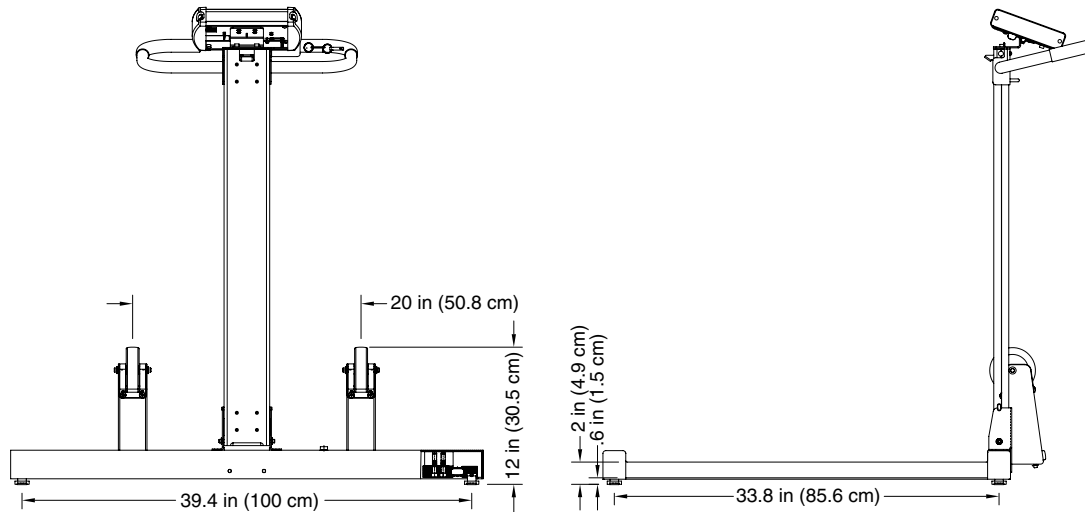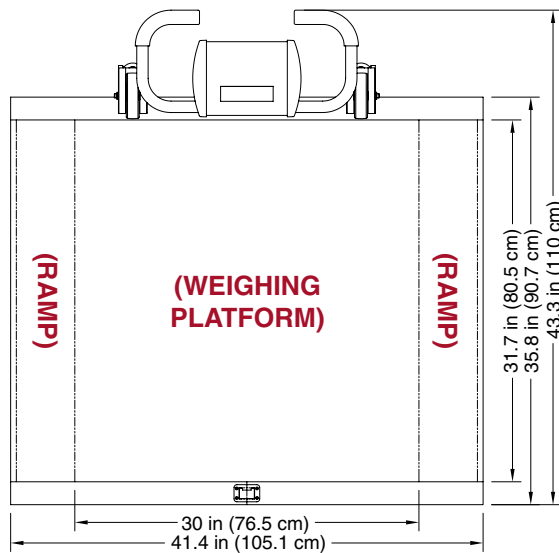

Capacity	1,000 lb x 0.2 lb / 450 kg x 0.1 kg
Display	Weight: Five digit, seven segment .5 in / 13 mm H LCD Height: Four digit, seven segment .375 in / 10 mm H LCD
Keyboard	Membrane type with 19 keys
Connectivity Ports	RS232 serial, USB, and optional Ethernet
Display Dimensions	9.5 W x 6.6 H x 2.1 D in / 24 W x 17 H x 5 D cm
Power	6 C cell Ni-Cad or NiMH batteries (not included) OR an optional 100 to 240 VAC 50/60Hz 12 VDC 1A wall plug-in UL/CSA listed AC power adapter (Detecto part number 6800-1045)
Resolution	5,000 divisions, maximum
Operating Temperature	14° to 104° F (-10° to +40° C)
Enclosure Construction	Steel and Aluminum
UPC Code	809161139402
Shipping Weight	91 lb / 41.2 kg

Detecto Scale reserves the right to improve, enhance, or modify features and specifications without prior notice.

 **DETECTO**[®]
A Division of Cardinal Scale Manufacturing Co.
203 E. Daugherty, Webb City, MO 64870 USA
Ph: 417-673-4631 or 1-800-641-2008 Fax: 417-673-5001
www.Detecto.com

Sold By:

