

MODEL EU544 DIGITAL LOW-PROFILE SCALE



- 200 kg x 0.1 kg capacity
- 12.7 mm high LCD
- compact design

**RECHARGEABLE
BATTERY POWER!**

Detecto's low-profile EU544 digital scale provides everything from basic patient weighing to printing or sending patient ID and weight data to a PC. With Body Mass Index calculation, automatic display lock, pushbutton and keypad tare, the EU544 is excellent for doctors' offices, hospitals, weight loss clinics, and long-term care facilities.



The New EU544


Low-Profile Digital Scale



SERIAL CONNECTIVITY PROVIDES:
Output To Optional Printers
Send Data To A PC

SIMPLE TO USE

For weighing only use, the yellow "quick keys" (pictured to the right) provide fast, easy operation, even by first time users. The full complement of displayed information makes it easy to read the data at a glance and the display hold function can be set for manual operation or automatic activation on stable weight.

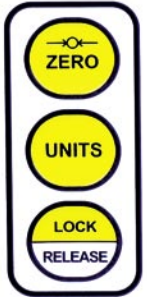

The new serial port allows output of the EU544's data directly to a printer by simply pressing the PRINT/ENTER key.



PATIENT ID

Enter up to an 11 digit patient ID that can be printed on a label or ticket or sent with the weight data to a PC. Detecto can provide the optional P150 battery-operated mobile ticket printer (pictured above) to professionally document your patients' weight by time and date for medical files. You can also print the patient height if it is required for their records.



The EU544 uses rechargeable Nickel Metal Hydride batteries, so you can take the scale to the patient. Optional AC adapter available.

BODY MASS INDEX

Body Mass Index (BMI) measurements have been shown to be an accurate estimate of health risk when used as part of a comprehensive



health assessment. The EU544 quickly calculates BMI by simply entering the patient's height while weighing.



Cardinal Scale reserves the right to improve, enhance or modify features and specifications without prior notice. All registered trademarks are the property of their respective owners.

EU544 SPECIFICATIONS	
Display:	Liquid Crystal, Reflective Type
Weight:	5 digit, 7 segment 13 mm high
BMI:	3 digit, 7 segment 13 mm high
Height-ID:	11 digits entered - 4 digits displayed, 7 segment 8 mm high
Bar Graph:	8 segment 8 mm high
Enclosure:	Stainless steel
Keypad:	19 key, color-coded, sealed membrane
BMI Range:	15.0 - 60.0
BMI Height Range:	90 cm to 215 cm
Capacity:	200 kg x 0.1 kg
Resolution:	5,000 divisions, maximum
Accuracy:	Class III
Power:	6 "C" size rechargeable NiMH cells (not included) or optional AC adapter.
Platform Size:	26 cm x 34 cm

SOLD BY:



DETECTO
A Division of Cardinal Scale Manufacturing Co.
203 E. Daugherty, Webb City, MO 64870 USA
Ph: 417-673-4631 or 1-800-641-2008 Fax: 417-673-5001
www.DetectoScale.com